



## ADVACRYL PLUS

### Polyglactin - 910

Antibacterial with Triclosan Braided Coated Absorbable Surgical Suture

**Consistent performance for confidence in a variety of soft tissue approximation minimizing Surgical Site Infections (SSI)**



**ADVA TRAY™ - the confidence to suture with ease**

ADVACRYL PLUS

Polyglactin - 910  
Antibacterial with Triclosan Braided Coated Absorbable  
Surgical Suture



With great versatility comes the confidence to meet specific needs

Product Description
















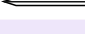

Sr. No.	Code	Suture			Needle				Units Per Box
		Colour	Size	Length (cm)	Curvature	Type	Length (mm)	Speciality	
1	ADVA247PH	VIOLET	0	70	1/2 CIRCLE		31	-	36
2	ADVA518PH	VIOLET	0	90	1/2 CIRCLE		36	-	36
3	ADVA946P	UNDYED	0	90	1/2 CIRCLE		36	-	12
4	ADVA2346P	VIOLET	0	90	1/2 CIRCLE		40	-	12
5	ADVA2517P	VIOLET	0	90	1/2 CIRCLE		36	-	12
6	ADVA2518P	VIOLET	0	90	1/2 CIRCLE		40	-	12
7	ADVA2534P	VIOLET	0	90	1/2 CIRCLE		36	-	12
8	ADVA519P	VIOLET	1	90	1/2 CIRCLE		36	-	12
9	ADVA947P	UNDYED	1	90	1/2 CIRCLE		36	-	12
10	ADVA2347P	VIOLET	1	90	1/2 CIRCLE		40	HEAVY	12
11	ADVA2347VMP	VIOLET	1	110	1/2 CIRCLE		40	HEAVY	12
12	ADVA2350P	VIOLET	1	120	1/2 CIRCLE		36	-	12
13	ADVA2359P	VIOLET	1	70	1/2 CIRCLE		40	HEAVY	12
14	ADVA2364P	VIOLET	1	90	1/2 CIRCLE		40	HEAVY	12
15	ADVA2421P	VIOLET	1	90	1/2 CIRCLE		40	HEAVY	12
16	ADVA2826P	VIOLET	1	35	1/2 CIRCLE		23	HEAVY	12
17	ADVA2826EP	VIOLET	1	70	1/2 CIRCLE		23	HEAVY	12
18	ADVA9335PH	VIOLET	1	70	J TYPE		31	-	36
19	ADVA9365PH	VIOLET	1	70	1/2 CIRCLE		36	-	36
20	ADVA2452P	VIOLET	2	90	1/2 CIRCLE		40	HEAVY	12
21	ADVAS2452P	VIOLET	2	90	1/2 CIRCLE		45	HEAVY	12
22	ADVA2478P	VIOLET	2	90	1/2 CIRCLE		40	HEAVY	12
23	ADVA2478PL	VIOLET	2	100	1/2 CIRCLE		40	HEAVY	12
24	ADVA9378PH	VIOLET	2	70	1/2 CIRCLE		45	-	36
25	ADVA945SP	UNDYED	2-0	90	1/2 CIRCLE		36	-	12







Round Bodied

Reverse Cutting

Cutting

Taper Cut

Sr. No.	Code	Suture			Needle				Units Per Box
		Colour	Size	Length (cm)	Curvature	Type	Length (mm)	Speciality	
26	ADVA2317P	VIOLET	2-0	90	1/2 CIRCLE		30	-	12
27	ADVA2382P	VIOLET	2-0	90	1/2 CIRCLE		40	-	12
28	ADVA2404P	UNDYED	2-0	90	3/8 CIRCLE		26	-	12
29	ADVA9375PH	VIOLET	2-0	70	1/2 CIRCLE		45	-	36
30	ADVA311PH	VIOLET	3-0	70	1/2 CIRCLE		22	-	36
31	ADVA316PH	VIOLET	3-0	70	1/2 CIRCLE		26	-	36
32	ADVA452PH	VIOLET	3-0	70	3/8 CIRCLE		24	-	36
33	ADVA458PH	UNDYED	3-0	70	1/2 CIRCLE		22	-	36
34	ADVA2236P	VIOLET	3-0	75	1/2 CIRCLE		36	-	12
35	ADVA2427P	VIOLET	3-0	90	1/2 CIRCLE		22	-	12
36	ADVA2437P	VIOLET	3-0	70	1/2 CIRCLE		20	-	12
37	ADVA397PH	VIOLET	4-0	70	3/8 CIRCLE		19	-	36
38	ADVA2304P	VIOLET	4-0	70	1/2 CIRCLE		20	-	12
39	ADVA2304PL	VIOLET	4-0	100	1/2 CIRCLE		20	-	12
40	ADVA2494P	VIOLET	4-0	45	3/8 CIRCLE		12	PRIME	12
41	ADVA9410PH	VIOLET	4-0	45	STRAIGHT		13	-	36
42	ADVA2303P	VIOLET	5-0	45	1/2 CIRCLE		16	-	12

Sr. No.	Double Arm Codes	Suture			Needles (1 & 2)				Units Per Box
		Colour	Size	Length (cm)	Curvature	Type 1	Type 2	Length (mm)	
43	ADVA2346DA180P	VIOLET	0	180	1/2 CIRCLE			40	12
44	ADVA2347DA140P	VIOLET	1	140	1/2 CIRCLE	 H		40	12
45	ADVA2347DA180P	VIOLET	1	180	1/2 CIRCLE	 H		40	12


 Round Bodied	 Reverse Cutting	 Cutting	 Taper Cut	 Oval Round Bodied	 Straight Cut
--	---	---	---	---	--

Product Description	
Composition	90% Glycolide , 10% L-lactide
Coating	Poly (gloycolide-co-lactide) (Glacomer 37) and Calcium Stearate, Triclosan
Size	U.S.P. Size 2 to 5-0
Colour	Undyed and Violet
Absorption Profile	Essentially complete up to 80 days
Needle Availability	A wide variety of Round Bodied (Taper Point), Ball Point, Cutting, Trocar, Reverse Cutting, Taper Cut, Blunt Point, Oval Round Bodied, Diamond Point & Spatula Micro Point needles

Tensile Strength	
Days of Implantation	Tensile Strength Retention
14 Days (2 To 4-0)	75%
21 Days (0 & 2-0 To 4-0)	55%
21 Days (2 & 1)	40%
28 Days (2 To 4-0)	15%

48% more knot pull tensile strength Vs. recommended limit ensures excellent support throughout the critical wound healing period

**Knot pull tensile strength (KgF)<sup>1</sup>**

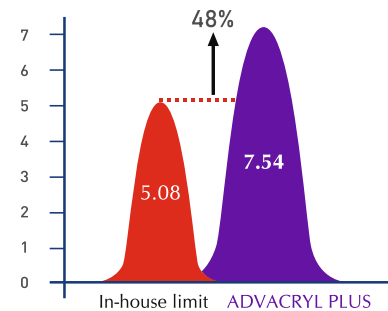

REFERENCE \_\_\_\_\_  
 1 : Bench test done on Triclosan coated Polyglactin-910 U.S.P. size -1, 40 mm Round bodied needle - Data on file with AMS.

### Multifilament Material



Provides excellent support for a range of critical wound healing period to secure closure in a variety of tissues, procedures and patients.



# ADVA TRAY™ - the confidence to suture with ease



## Premium Needles for consistent Superior Performance

**PREMIUM NEEDLES**  
ensure  
**UNSURPASSED SUTURING EXPERIENCE**

<b>DURABLE COAT</b> RETAINS SHARPNESS... PASS AFTER PASS	Maintains needle sharpness over the duration of procedure, resulting into consistent penetration pass after pass thus less potential for tissue trauma.
<b>RIBBED BODY</b> FOR SUPERIOR CONTROL	For superior grasping and stability of needle in needle forceps.
<b>CHAMFERING</b> MINIMUM SUTURE DETACHMENT	Diminishes chances for suture fracturing at the swage by rounding off the sharp edges of the barrel in the swage area of the needle.
<b>ELECTROPOLISHING</b> SMOOTH & BRIGHT	Makes needle surface free from micro asperities, thus creates ultra sharp edges for lesser tissue drag.

Minimal tissue trauma enables faster recovery